
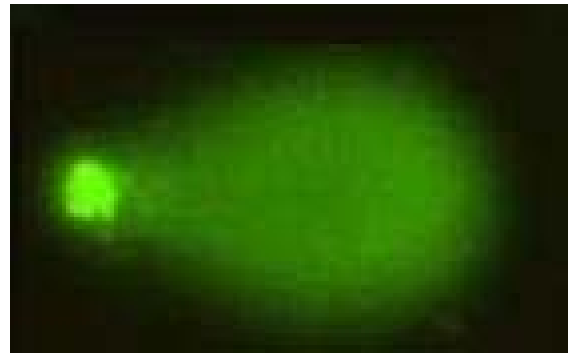


# TREVIGEN®

## *Standardize your comet assay....* CometAssay™ Electrophoresis System

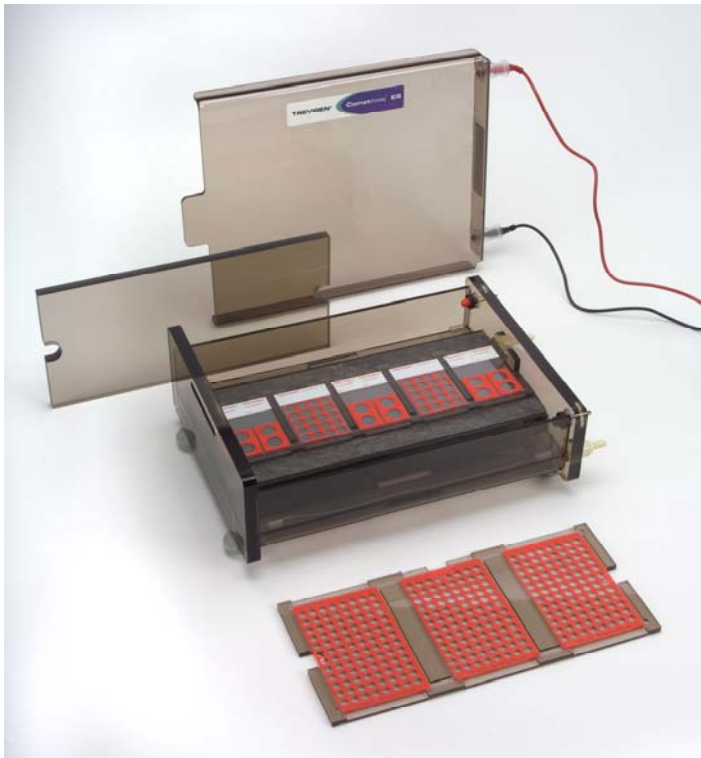

Trevigen's CometAssay™ ES overcomes variations by placing an acrylic overlay on top of an elevated slide tray to maintain optimal buffer height for DNA migration. A constant buffer temperature is maintained using an underlying water jacket to cool the ceramic slide platform and buffer chamber. Notice the outlet valves for the water jacket on the right side of image. Specially designed slide trays are provided to accommodate 2, 20 and 96 well slides and maintain proper slide orientation in an electrophoretic field uniquely configured for single cell gel electrophoresis.

**Cat#: 4250-050-ES**

***Coming Soon!***

### Features:

- Maintains constant buffer temperature
- Maintains optimal buffer level for consistent results
- Specially designed trays accommodate 2, 20 and 96 well slides and maintain correct position during electrophoresis
- Optimized for use with CometAssay™ Kits and CometAssay™ Control Cells
- Greater reproducibility and reduced variability between individual users and different labs

1-800-TREVIGEN (1-800-873-8443) • 301-216-2800 • Fax: 301-560-4973

Email: [sales@trevigen.com](mailto:sales@trevigen.com) • <http://www.trevigen.com>